



Vertical Needle Holder

Positions needles for easy loading

Conveniently positions 20 needles vertically in numbered positions for easy loading.


Loaded needles are shielded with 1/8" of stainless steel.

The tips of the needles are in no danger of blunting as needles are suspended. The needles are never resting on the fragile tip.


Can be used with any existing L-block shield.

Additional Needle Holders can be used to store loaded needles.

Needles are conveniently positioned for an autoradiograph.

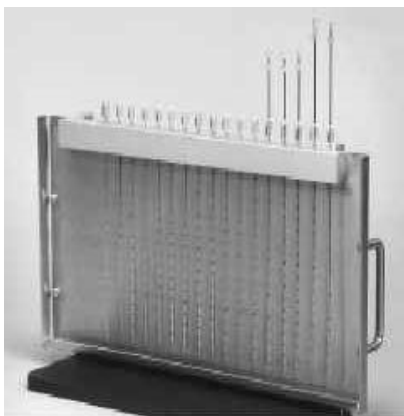

Needles in the Vertical Needle Holder ready for loading.

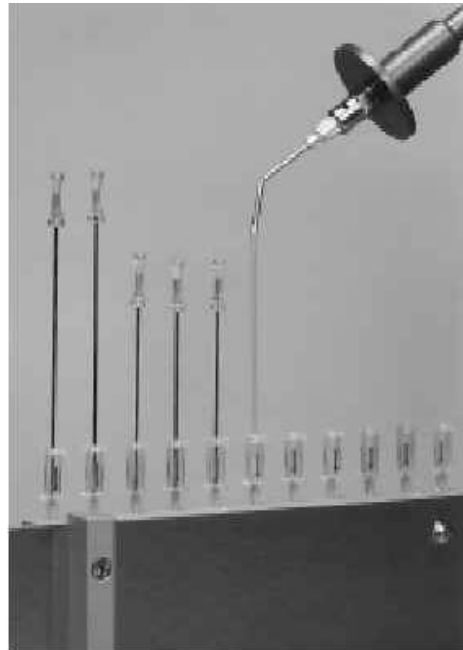

The SeedVac™ loading needles in the Vertical Needle Holder

© Elimpex-Medizintechnik, Spechtgasse 32, A-2340 Moedling, Austria
phone +43-2236-410450
fax +43-2236-410459

