




## SEED ALIGNMENT TRAY

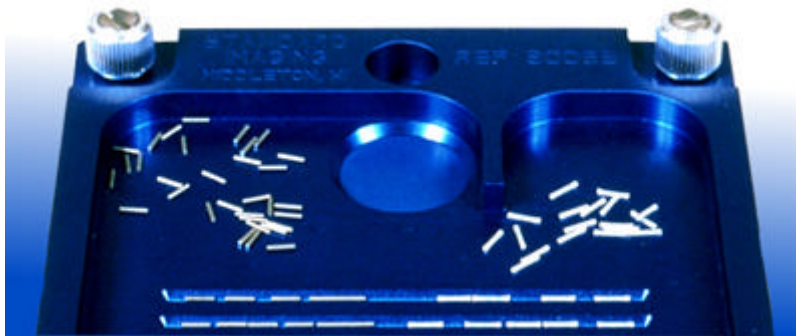

**Aligns seeds and  
spacers for easy access**

- Seeds can all be aligned in one direction for easy pick up by SeedVac™.
- Uniform alignment reduces arm movement.
- Uniform alignment speeds the entire needle loading procedure.
- Six grooves are provided for easy access of seeds or spacers.
- Interfaces with the Seed Sterilization and Sorting Tray for protection from radiation.
- All seeds are stored inside the tray and accessed as needed.
- Seeds are simply scooped out of the tray wells and slid into position using the provided spatula.



The Seed Alignment Tray connected to the Seed Sterilization and Sorting Tray. A SeedVac™ is picking up the seeds.



The Seed Alignment Tray with the seeds and spacers ready for loading.

